




DELAWARE EMERGENCY MEDICAL REPORTING SYSTEM (DEMRS)

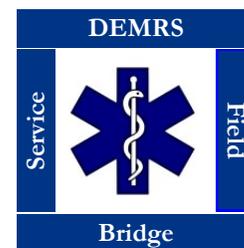

Installing ImageTrends EMS Field Bridge

Field Bridge Specifications and Licenses

Each agency will be issued one EMS Field Bridge License and Service Support Agreement from the Office of Emergency Medical Services (OEMS). The license key provided may be used to install Field Bridge on your agency's multiple laptops. The following specifications are necessary for Field Bridge to run optimally:

Required:	Operating Systems Supported:	Additional Software:	Optimal:
1GHz Processor (non-Atom)	Microsoft Windows XP, Vista, Windows 7, Windows 8	Microsoft .NET Framework 4.0 or higher	Multi-core processor; Intel Core i5 or greater
2 GB RAM		Adobe Flash Player 11 or higher	4 GB RAM or greater
600 MB Available Hard Disk Space	Screen Resolution Supported:	Adobe Reader 10 or higher	64-bit Windows operating system
802.11 Wireless	1024 x 768		1GB Available Hard Disk Space
			802.11 Wireless
			Bluetooth
			GPRS/CDMA Cellular modem

Downloading the Field Bridge Application

1. **Open** the web browser and **navigate** to www.imagetrend.com/support.
2. **Select** "EDS Application Downloads" from the menu.
3. A login page will display. **Enter** the username "DelawareOEMS" and Password "fieldbridge". **Select** "Login".

Login

Username:

Password:



4. **Select** the menu item "Field Bridge Version 5 Downloads". The EMS Field Bridge Version 5.x Downloads page will appear.

Installing ImageTrends EMS Field Bridge

5. **Scroll down** until you see Field Bridge v5 Full Install - Executable Version. **Read** through the notes and ensure prerequisites are met

Field Bridge v5 Full Install - Executable Version

The following is for complete installation of EMS Field Bridge on a machine that does not already have Field Bridge installed or requires a full re-installation.

Prerequisites: For users running Windows Vista or prior versions, install the .NET Framework 3.5 SP1 first, and then run the FieldBridgeInstaller.msi file.

ATTENTION: Do NOT install Field Bridge Version 5.2 or above unless your Service, State or Rescue Bridge is on version 5.2 or above.

Download Details	
Prerequisite:	.NET Framework 3.5 SP1
File Name:	setup.exe (489 Kb - Click green link below)
Additional Components:	Adobe Reader (used to view PDF reports) Adobe Flash Player (used to display the human assessment and power tools) PDF Creator (used for ZOLL and LIFEPAK EKG imports)
Date Published:	April 9, 2012
Version:	5.2.5.2 Setup Instructions


[Download Executable Field Bridge 5 Installer](#)

6. **Click** on the green arrow “Download Executable Field Bridge 5 Installer”


7. A setup.exe file prompt will launch and a File Download - Security Warning, **select Run.**




8. A secondary Internet Explorer - Security Warning prompt may display, **select Run.**



9. Allow the Image Trend EMS Field Bridge 5 setup to complete, then the Setup Wizard will display, **select Next**

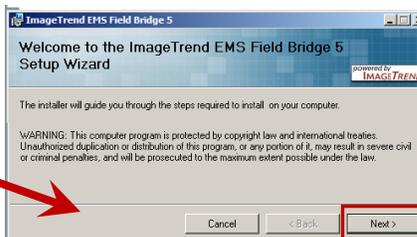

10. You will be prompted to Select Installation Folder. Select a folder or allow the application to install to the local drive. **Select Next** to continue the installation process.

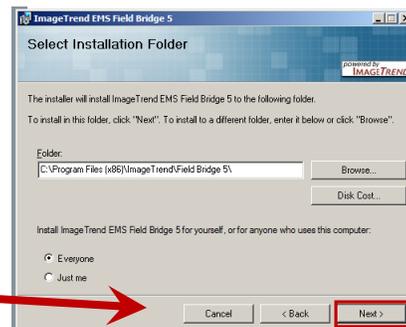

Installing ImageTrends EMS Field Bridge

11. Once installation is complete, **select** Close to exit. You can now close your browser window.





Field Bridge 5 icon will appear on the desktop. You have completed Field Bridge 5 Installation.






Registering the Field Bridge Application

1. **Select** the Field Bridge 5 Icon to launch the application. The first time the application launches you will need to Enter Registration Key. The registration key is a 25 digit code found on EMS Field Bridge License and Service Support Agreement that was issued to your Agency by OEMS.
2. **Enter** the 25 digit code.
3. **Select** OK.
4. The EMS Field Bridge login prompt will appear. **Select** Settings.

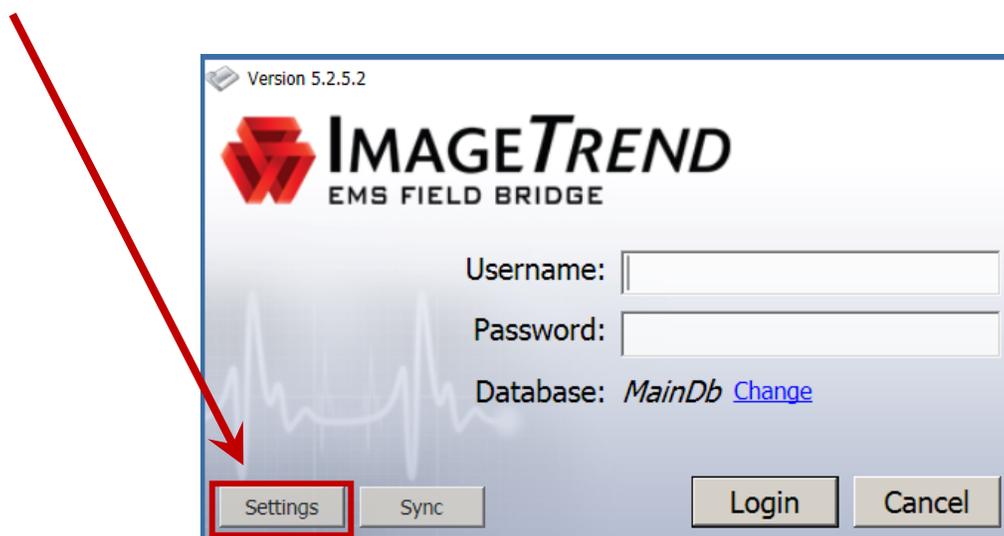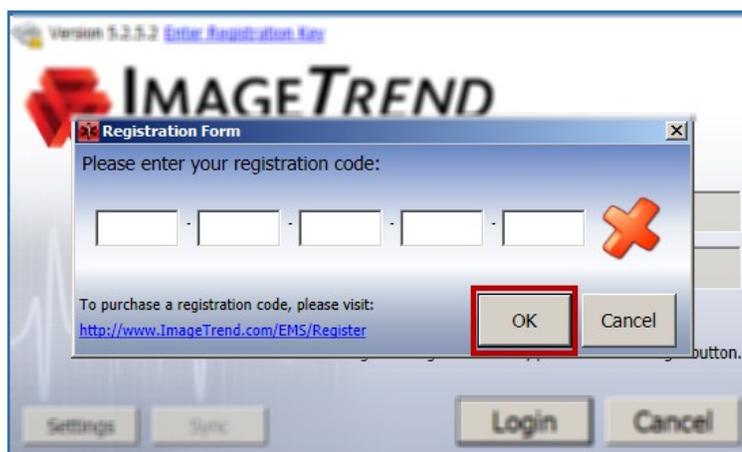


Installing ImageTrends EMS Field Bridge


5. On the Settings Sync/Post screen, in the Agency ID field **enter** your Station Number.
6. On the Settings Sync/Post screen, in the Service URL field, **enter** :
 https://delaware.emsbridge.com
7. **Select** OK.
8. In order to complete the setting changes, you will be prompted to Validate Sync Change. Use your personally assigned Field Bridge Username and Password, **select** OK.

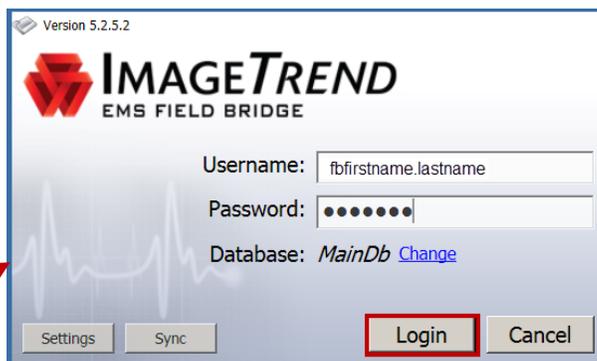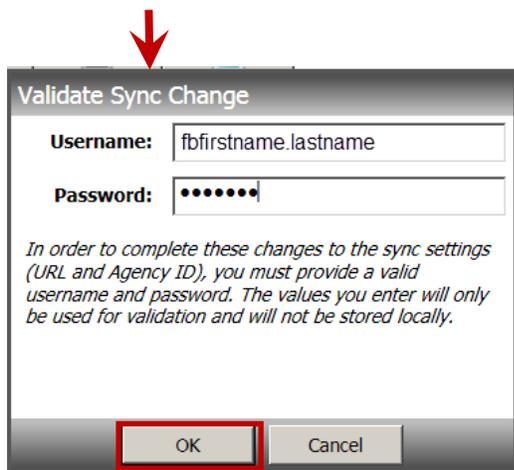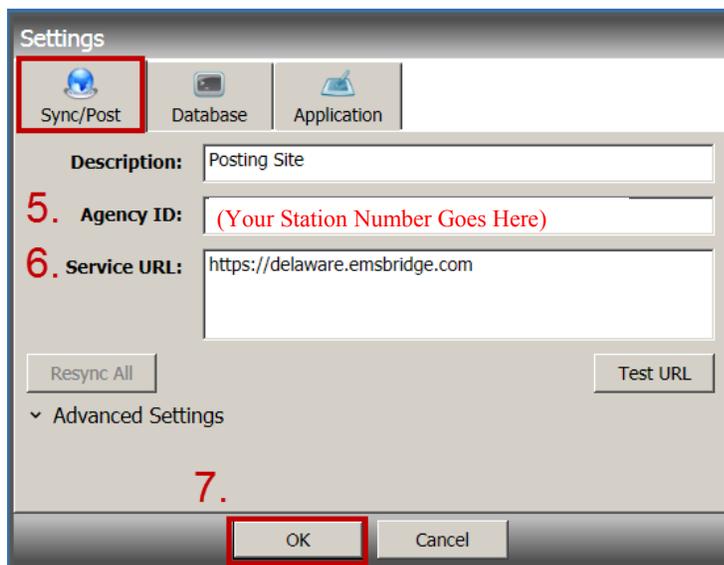

9. Once Field Bridge Sync completes, the Field Bridge login prompt will appear. Using your personally assigned Field Bridge Username and Password **Login** to Field Bridge.
10. The Field Bridge application will launch and notify you that Updates are Available. **Select** Yes.

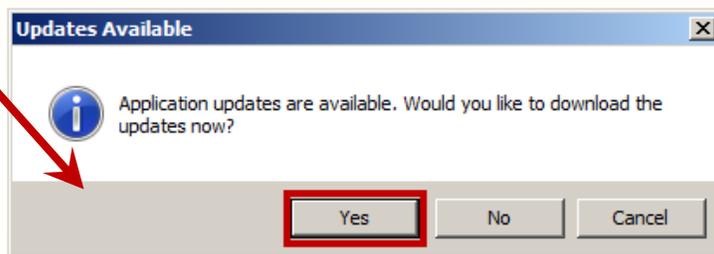


Installing ImageTrends EMS Field Bridge

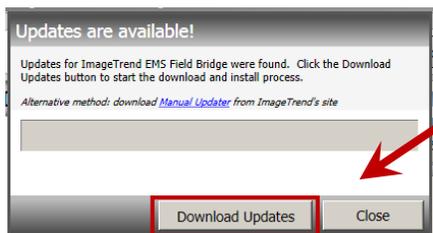


11. Again you will be prompted that Updates are available, **select** Download Updates.

12. The updates will download and the message box will automatically close.

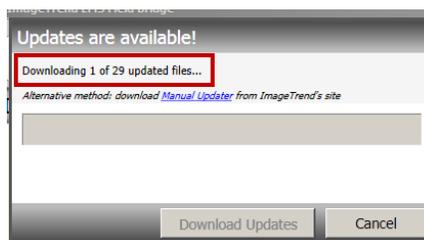

13. **Select** Log Out to close the Field Bridge Application.

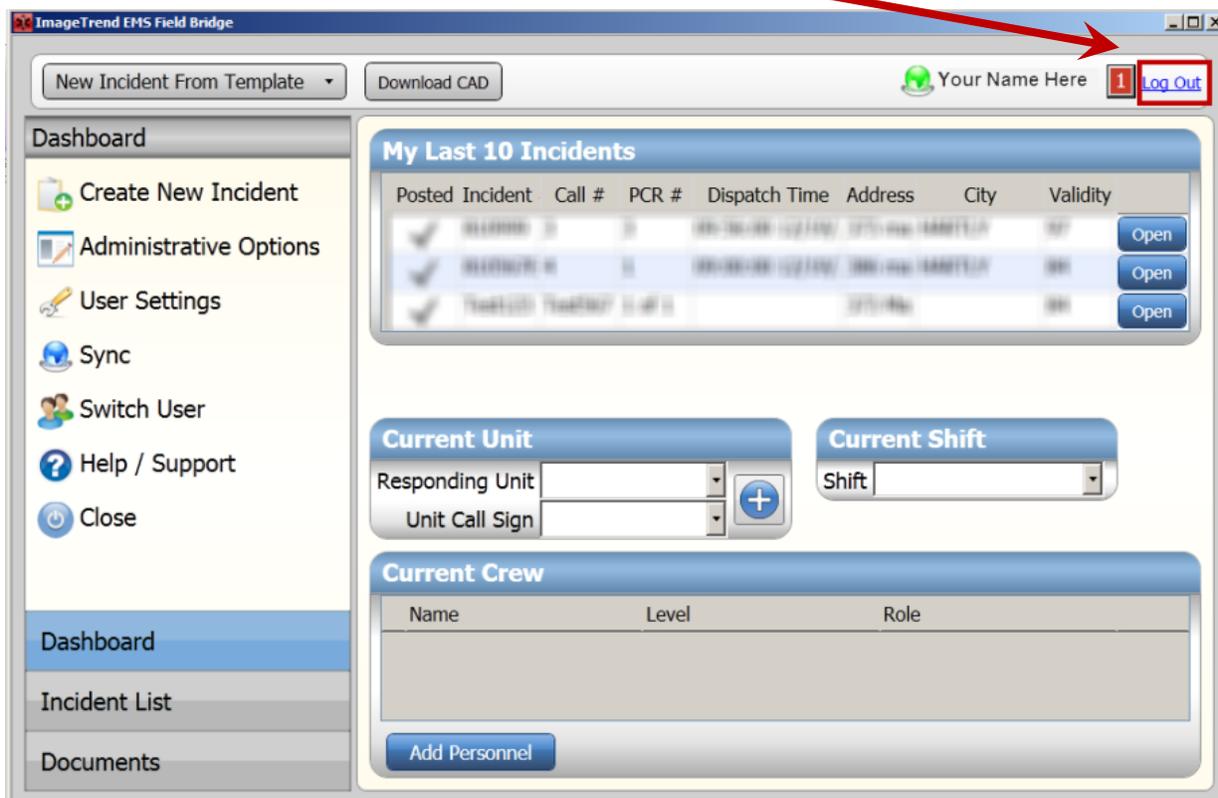

Registration and updates to Field Bridge 5 are complete.